



Product information **sc5010 SEK**



Mobile surgical chair



Here you will find the technical data
for the mobile surgical chair

sc5010 SEK



Mobile surgical chair

The mobile operating chair for large joints

The mobile operating chair sc5010 SEK (shoulder, elbow, knee) allows outpatient procedures on large joints to be carried out. These include operations in the areas of orthopedics and arthroscopy. The flexible structure of the sc5010 SEK operating chair also allows it to be used in the fields of dermatology and aesthetic plastic surgery.

Adapted to the representation of large joints, the height adjustment of the seat from 63cm to 98cm, as well as the individual adjustment options for the head, back, torso and leg supports, support the precise adjustment of the operating field. Standard lifting columns (up to 200kg patient weight) or heavy-duty lifting columns (up to 300kg patient weight) are available. Backrest segments can be easily removed if necessary, providing unrestricted access to the shoulder joint. The flexibly adjustable leg supports can be adjusted horizontally and vertically or removed completely. Various headrest systems are available for the required head support or fixation.

The clearly arranged control unit has an on/off switch, a battery capacity display and buttons for controlling the motorized seat, back and head segments. Six memory functions make it easy to quickly access repetitive settings. In an emergency, the lying surface can be brought into the shock position at the push of a button. Manufactured AKRUS cushions ensure a relaxed sitting and lying feeling during the treatment. A wide selection of colors and two upholstery thicknesses are available.

Matching accessories round off the functionality of the sc5010 SEK.

The mobile operating chair for large joints sc5010 SEK is delivered fully assembled and the use of high-quality components allows many years of reliable use when used as intended. In the event of a service requirement, authorized medical technology dealers are available.



Surgical table, flat

*Chair configuration to
Versatile leg adjustments*

Deck chair configuration


Shock position

Focus: Patient positioning and accessibility

APPLICATIONS



*Removable shoulder segments
Option: TCHR – Head fixation
headrest with inclinable holder*

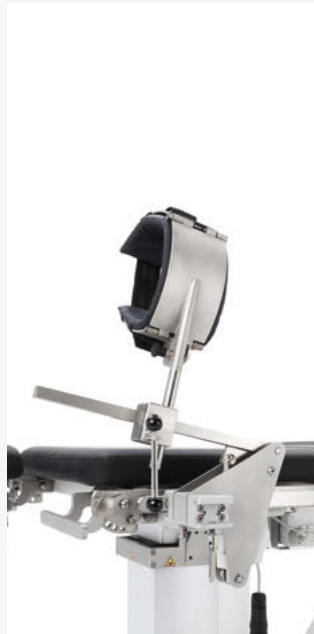

Knee holder systems



Goepel leg supports



Forearm extension table



Flexible patient positioning

Application areas



Shoulder



Knee



Hand / Arm



Foot / Leg

ACCESSORIES

*Forearm extension table**Manual knee holder**Electric knee holder**Goepel leg supports**Adjustable head board**Multi articulated head rest, inclinable**IV pole supports**Drape support, solid**4 Function foot joystick panel**Foot switch 2 axes**Battery**Control unit*

Upholstery fabrics

Choose one of 16 colors from our standard range. Special colors on request.




Combra White /
207x605




Creme /
208x2380




Sand /
207x3261

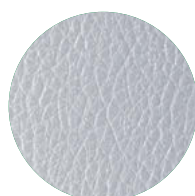


Pearl /
032x4354

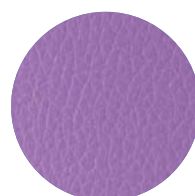

Lilac /
208x2392

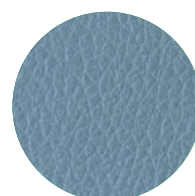


Bluegrey /
208x1277

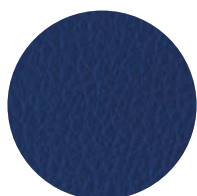


Corn flower /
207x729



Midnightblue /
032x4709

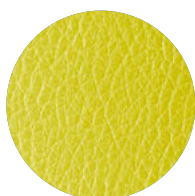

Yellow /
208x2795

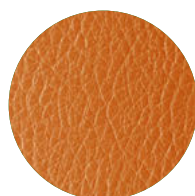


Orange /
208x3107




Red /
207x2796




Oper /
218x1637

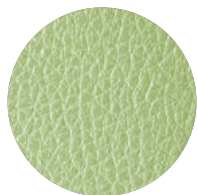


Mint /
208x2394

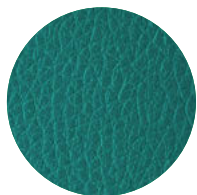

Turquoise /
208x749

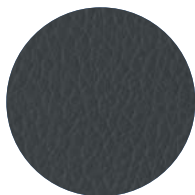

Charcoal /
207x1049

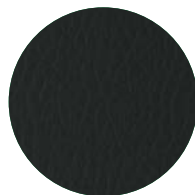

Black /
electrically conductive
anti-static
217x3117

Distributor