

Product information
ak480



The flexible surgeon's chair



Here you will find the technical data for the surgeon's chair

The flexible surgeon's chair

Electrically height-adjustable surgeon's chair

The electrically height-adjustable surgeon's chair ak480 is the solution for seated surgical procedures in the fields of ophthalmology, neurosurgery and head surgery. Thanks to its versatile adjustment options, the chair offers fatigue-free, comfortable sitting and maximum support for the surgeon during the procedure.

With five wheels and mains-independent battery operation, the ak480 surgeon's chair is very mobile. The height of the seat is adjusted easily using two separate heel pedals. Based on the surgeon's preferences, the ergonomically shaped seat can be tilted and moved horizontally, depending on the equipment. A backrest with height and tilt adjustments as well as lumbar support allows fatigue-free sitting. The AKRUS chair fixation system ensures stability during the procedure. The chair is blocked by pressing on the easy-to-reach foot ring. To unlock, press the middle foot plate with your heel.

The armrest is also designed to be flexible and fatigue-free. The pre-shaped multi-directional multi-joint armrest offers unlimited adjustment options to adjust any ergonomically desired or operationally required position. They are easy to navigate with your elbows to adjust the optimal position or get out of the chair.

The ak480 surgeon's chair is delivered fully assembled and the use of high-quality components allows many years of reliable use when used as intended. In the event of a service requirement, authorized medical technology dealers are available.



Swiveling arm rests



Moveable armrest holder



Seat and backrest tilt



Anatomical seat cushion

Focus: ergonomic comfortable seating

APPLICATIONS



Armrests swivel: in



Armrest swivel: out



Seat tilt and backrest tilt



Heel controls: seat up, roller fixing, seat down

Application areas



Ophthalmology



ENT, OMF



Plastic, aesthetic
and neurosurgery



Dental

FUNCTIONS



Multi articulated armrests shifting: fwd



Multi articulated armrests shifting: bwd



Multi articulated armrests tilting



Seat shift: fwd



Seat shift: bwd



Seat tilt: down



Seat tilt: up



Heel pedals



Multi articulated armrest



Anatomical seat cushion



Backrest with lumbar support

The ak480 surgeon's chair can be ordered with standard armrests, standard backrest and fixed seat holder. Mind the reduced functionality when placing the order.

Upholstery fabrics

Choose one of 16 colors from our standard range. Special colors on request.




Combra White /
207x605




Creme /
208x2380




Sand /
207x3261

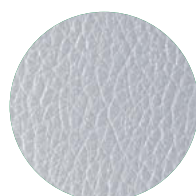


Pearl /
032x4354

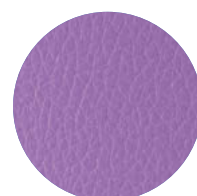

Lilac /
208x2392

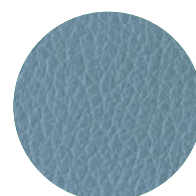


Bluegrey /
208x1277

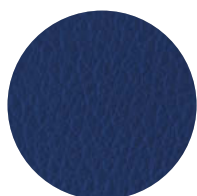

Corn flower /
207x729



Midnightblue /
032x4709

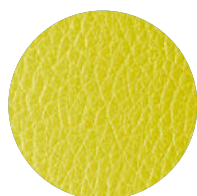

Yellow /
208x2795




Orange /
208x3107




Red /
207x2796




Oper /
218x1637

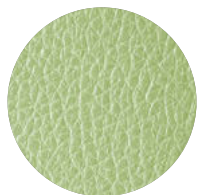


Mint /
208x2394

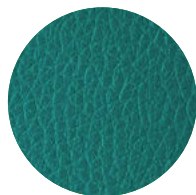

Turquoise /
208x749

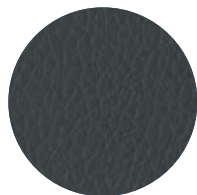

Charcoal /
207x1049

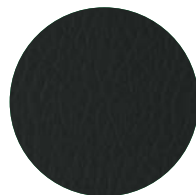

Black /
electrically conductive
anti-static
217x3117

Distributor