

Product information  
**sb5010 ES/HS**



## The mobile operating table



Here you will find the technical data  
for the mobile surgical chairs

sb5010 ES



sb5010 HS



## Mobile surgical bed

### The mobile operating table for versatile outpatient applications

The mobile, battery-operated operating table series sb5010 ES/HS was developed for versatile outpatient applications, including: ophthalmology, ENT (ear, nose and throat), OMF (oral and maxillofacial) as well as plastic, reconstructive and aesthetic surgery. Compared to the mobile operating chairs, this table has a surface area that is 10cm longer.

Adapted to outpatient treatments in the head area, the height adjustment of the mobile operating table **sb5010 HS (Head Surgery)** from 61cm to 96cm, as well as the individual adjustment options for the head, back and torso elements, support the precise adjustment of the operating field. Standard lifting columns (up to 200kg patient weight) or heavy-duty lifting columns (up to 300kg patient weight) are available. Quality rollers and the AKRUS roller fixation system ensure mobility and a secure stand.

The special need for a low ingress level and comfortable leg room for the operating surgeon in the field of eye surgery is answered by the mobile operating table **sb5010 ES (Eye Surgery)**. Its low ingress level of just 52cm allows specifically elderly patients to enter the bed with a great feeling of safety and comfort. The bed lifts up to 78cm and allows precise adjustment options for the operating field.

The clearly arranged control unit has an on/off switch, a battery capacity display and buttons for controlling the motorized lounger, back and headrest segments. Six memory functions make it easy to quickly access repetitive settings. In an emergency, the bed surface can be brought into the shock position at the push of a button. Manufactured AKRUS cushions ensure a relaxed sitting and lying feeling during the treatment. A wide selection of colors and two upholstery thicknesses are available. Matching accessories round off the functionality. AKRUS mobile operating tables are delivered fully assembled and the use of high-quality components allows for many years of reliable use when used as intended. In the event of a service requirement, authorized medical technology dealers are available.



*Very low ingress level*



*Deck chair configuration*



*Spacious leg room under backrest*



*Shock position*

## Focus: Head positioning

### APPLICATION PROPOSALS



*Ophthalmology*



*ENT surgery*



*OMF surgery*



*Individual patient positioning*

## Application areas



Ophthalmology



ENT, OMF



Plastic and  
aesthetic surgery

## ACCESSORIES



*Std. head rest incl. surgeon wrist rest*



*Multi articulated head rest, tiltable holder*




*Padded head board, tiltable holder*



*Safety railing*




*Arm rest, multi articulated*

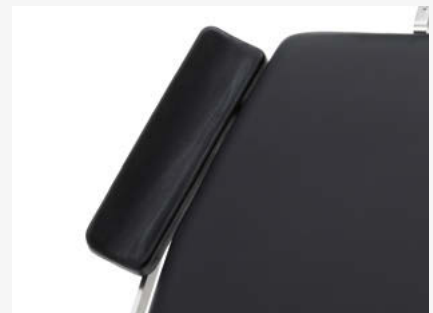

*Arm rest, standard*



*IV pole*

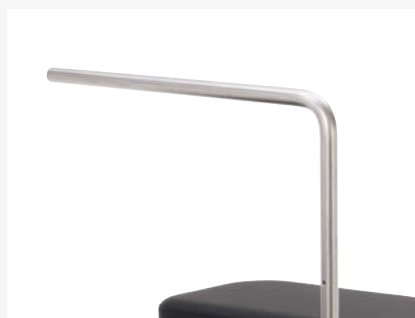

*Drape support, solid*



*Control unit*



*Foot controls for castors*



*Foot switches*



*Foot Joystick panel*

## Upholstery fabrics

*Choose one of 16 colors from our standard range. Special colors on request.*



Combra White /  
207x605




Creme /  
208x2380




Sand /  
207x3261




Pearl /  
032x4354

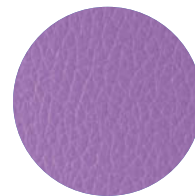

Lilac /  
208x2392

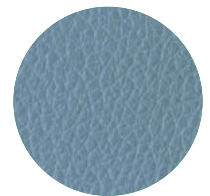


Bluegrey /  
208x1277

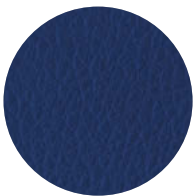

Corn flower /  
207x729



Midnightblue /  
032x4709

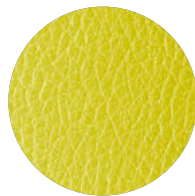

Yellow /  
208x2795




Orange /  
208x3107




Red /  
207x2796




Oper /  
218x1637

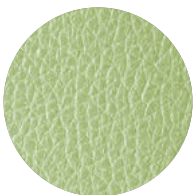


Mint /  
208x2394

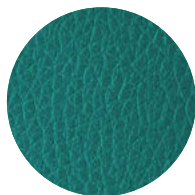

Turquoise /  
208x749

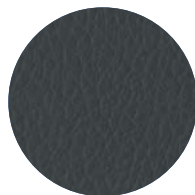

Charcoal /  
207x1049

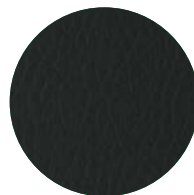

Black /  
electrically conductive  
anti-static  
217x3117

Distributor